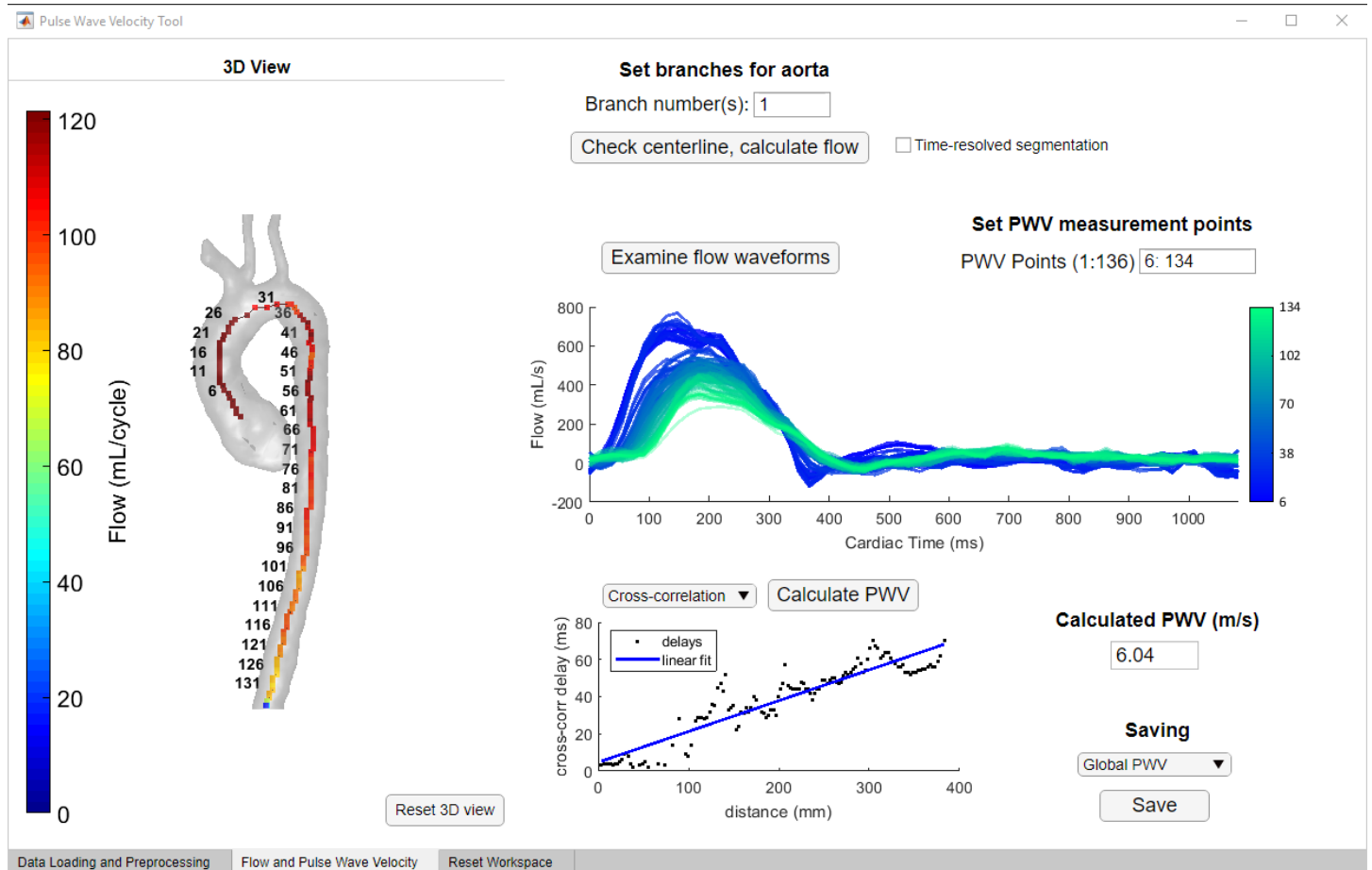


## 4D Flow Pulse Wave Velocity Tool – MATLAB-based App

Date: March 31, 2021

Author: Eric M. Schrauben, PhD

[e.m.schrauben@amsterdamumc.nl](mailto:e.m.schrauben@amsterdamumc.nl)

This software builds from previous work, and the use of this tool should be associated and cited with the following reference:

E. Schrauben, et al. "Fast 4D flow MRI intracranial segmentation and quantification in tortuous arteries." JMRI 2015, doi: 10.1002/jmri.24900.

### Installation:

This app is built using MATLAB's appdesigner functionality. It was built using Matlab version 2019a; it is therefore recommended to use this version or newer for running this app. Download and unzip all files. All functionality of the app is built into PulseWaveVelocityTool.mlapp and the corresponding utilities subfolder.

### Data needed

The app works directly on reconstructed .par / .rec files from a Philips 4D flow acquisition. It assumes these are all located within a single folder. The files needed are as follows:

- SUBJECTID or other identifier\_1.par, SUBJECTID or other identifier\_1.rec
- SUBJECTID or other identifier\_2.par, SUBJECTID or other identifier\_2.rec
- SUBJECTID or other identifier\_3.par, SUBJECTID or other identifier\_3.rec

These correspond to the magnitude and phase images for each velocity direction (vx, vy, vz).


## Basic Processing Steps and Information

### 1. Getting Started

Within Matlab, navigate to the folder containing PulseWaveVelocityTool.mlapp. Type this into your command window:

```
>> PulseWaveVelocityTool
```

Here is the initial screen:

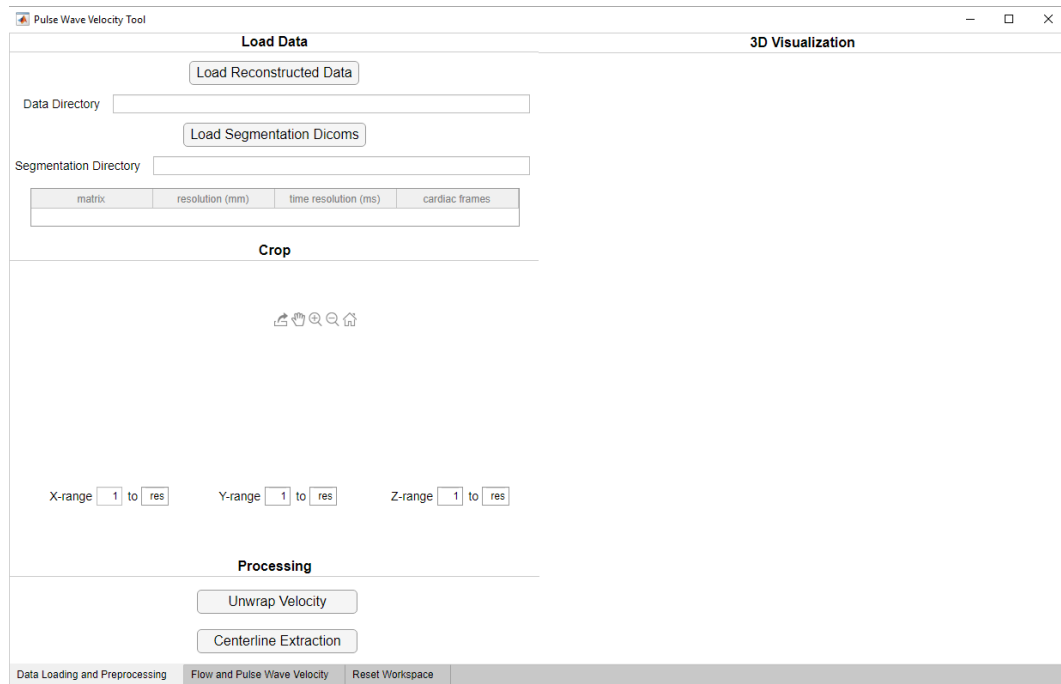


Processing is sorted into two main tabs: Data Loading and Preprocessing, and Flow and Pulse Wave Velocity


## 2. Loading and Preprocessing

Click Load Reconstructed Data. You will be prompted to find the one of the reconstructed .rec files in order to locate the correct data folder. Example aortic 4D flow data are provided for initial testing.


Following loading of data, a background subtraction GUI is automatically opened. Adjust the complex difference (CD) and Noise sliders to create a mask that only includes static tissue. Then click Update Images to calculate the background phase fit, and click Done to apply this fit to the velocity data.

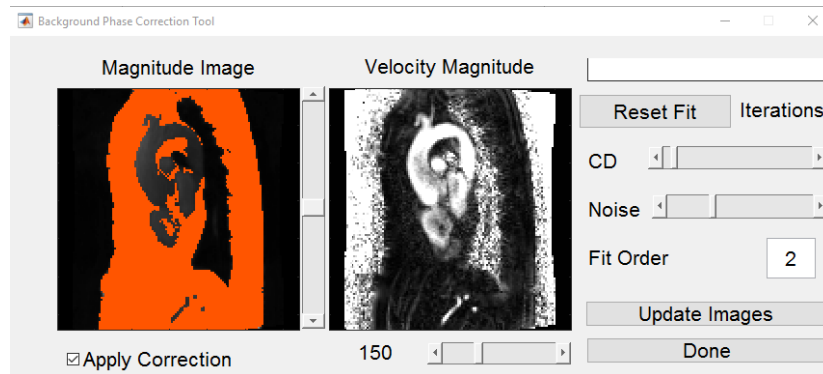

After loading completion, the tool will automatically display axial, coronal, and sagittal MIPs using a calculated angiographic image. It also automatically thresholds this image to produce a 3D isosurface within the field of view:

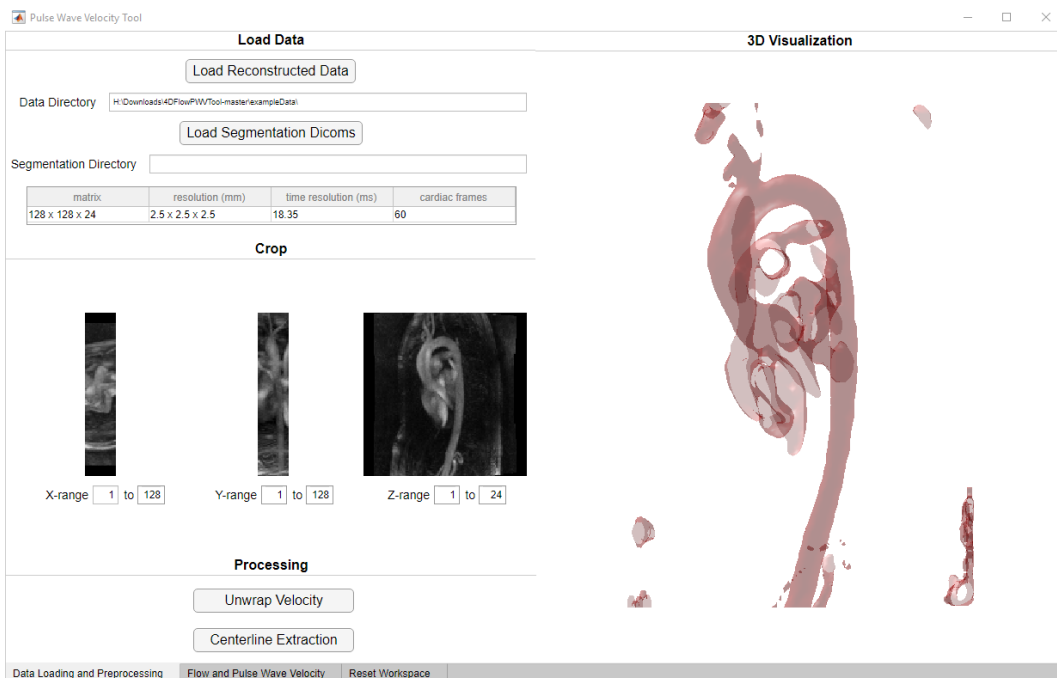

For best performance, a manual segmentation of the aorta should also be included. This segmentation is read in as dicoms produced from a given segmentation software. Click Load Segmentation Dicom to load and view the manual segmentation in gray overlaid on the previously calculated 3D isosurface:

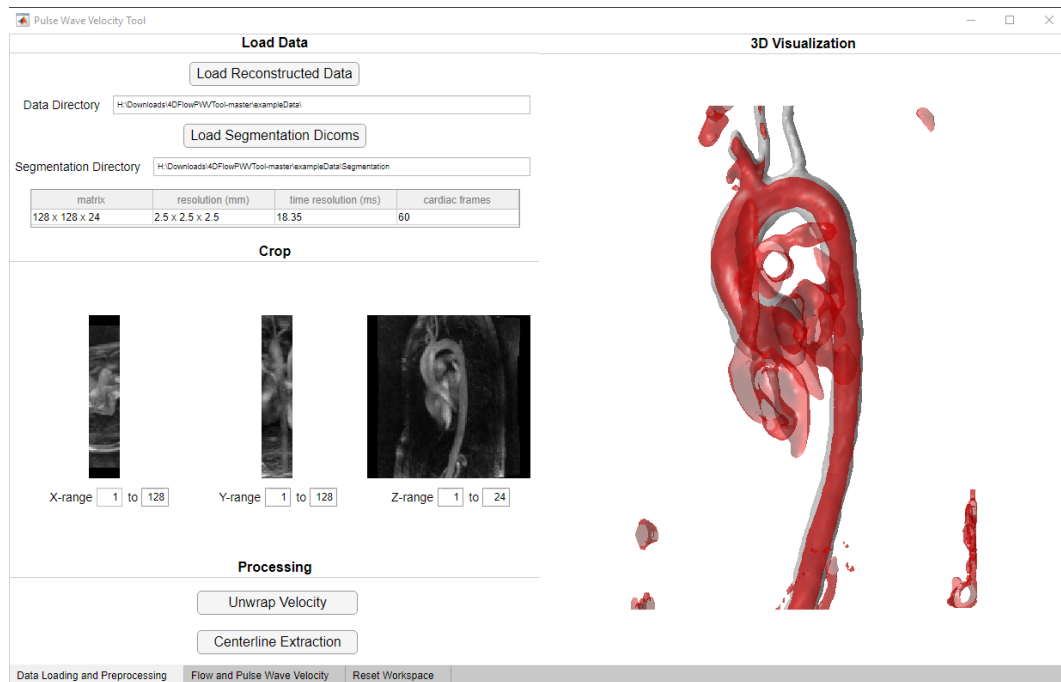

*Optional:* Automatic velocity unwrapping can be performed using the Unwrap Velocity button. For details please see and cite this work:


[J Magn Reson Imaging](#). 2016 Apr;43(4):833-42. doi: 10.1002/jmri.25045. Epub 2015 Sep 28.

### Phase unwrapping in 4D MR flow with a 4D single-step laplacian algorithm

Michael Loecher<sup>1</sup>, Eric Schrauben<sup>1</sup>, Kevin M Johnson<sup>1</sup>, Oliver Wieben<sup>1 2</sup>

Finally, to extract specific points within the aorta to perform pulse wave velocity calculation, push Centerline Extraction. The 3D vasculature is skeletonized and individual vessel centerlines are extracted and labelled with unique identifiers. Upon completion, the Flow and Pulse Wave Velocity tab is opened.

### 3. Flow and Pulse Wave Velocity

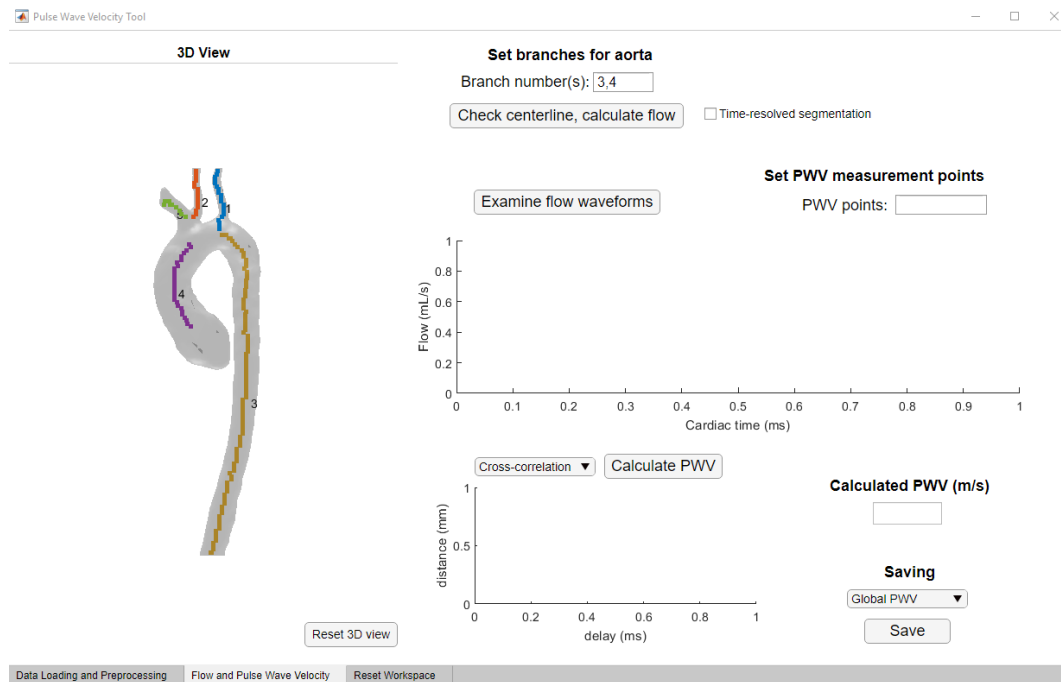

*Optional:* Time-resolved segmentation of the aorta can be performed by checking the Time-resolved segmentation checkbox. For details please see and cite this work:

[J Cardiovasc Magn Reson. 2017; 19: 47.](#)

Published online 2017 Jun 23. doi: [10.1186/s12968-017-0360-8](#)

PMCID: PMC5481950


PMID: [28645326](#)

**Improving visualization of 4D flow cardiovascular magnetic resonance with four-dimensional angiographic data: generation of a 4D phase-contrast magnetic resonance CardioAngiography (4D PC-MRCA)**

[Mariana Bustamante](#)<sup>1,2</sup> [Vikas Gupta](#)<sup>1,2</sup> [Carl-Johan Carlhäll](#)<sup>1,3</sup> and [Tino Ebbers](#)<sup>1,2</sup>

To segment the whole aorta, individual branches need to be selected. Input the corresponding branch numbers in Branch number(s):. Then click Check centerline, calculate flow. At each centerline point, an orthogonal cross section is calculated (using the norm of the centerline). Flow is calculated within each cross over all cardiac time frames. This step can take some time, depending on number of cardiac frames and the number of final centerline points within the aorta.

Following flow calculation, the total flow at each point and corresponding waveforms are plotted for inspection:

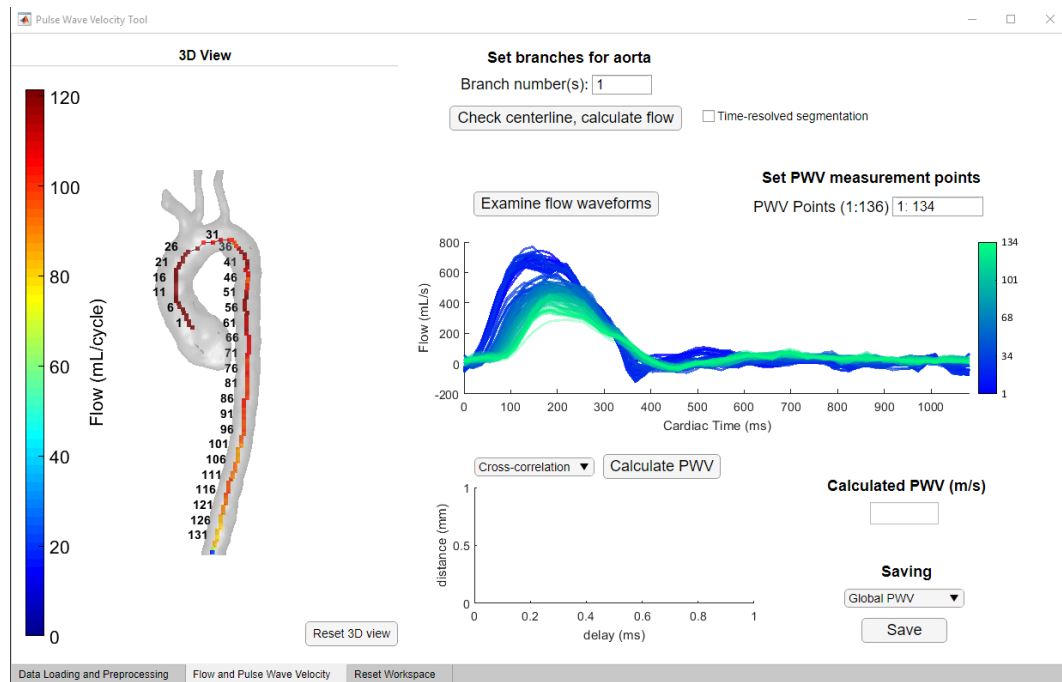

Finally, pulse wave velocity can be calculated over any number of points along the aorta. Input these points into PWV Points, then click Calculate PWV. The calculation of PWV has three options:

- Cross-correlation:

[J Magn Reson Imaging. 2008 Jun;27\(6\):1382-7. doi: 10.1002/jmri.21387.](#)

### A new method for the determination of aortic pulse wave velocity using cross-correlation on 2D PCMR velocity data

Samuel W Fielden<sup>1</sup>, Brandon K Fornwalt, Michael Jerosch-Herold, Robert L Eisner, Arthur E Stillman, John N Oshinski

- Wavelet cross spectrum analysis:

[Comparative Study](#) [J Cardiovasc Magn Reson. 2015 Jul 30;17\(1\):65.](#)


doi: 10.1186/s12968-015-0164-7.

### Estimation of aortic pulse wave transit time in cardiovascular magnetic resonance using complex wavelet cross-spectrum analysis

Ioannis Bargiotas<sup>1</sup>, Elie Mousseaux<sup>2,3</sup>, Wen-Chung Yu<sup>4</sup>, Bharath Ambale Venkatesh<sup>5</sup>, Emilie Bollache<sup>6</sup>, Alain de Cesare<sup>7</sup>, Joao A C Lima<sup>8</sup>, Alban Redheuil<sup>9,10,11</sup>, Nadja Kachenoura<sup>12</sup>

- Time-to-foot

The output of this calculation is shown:

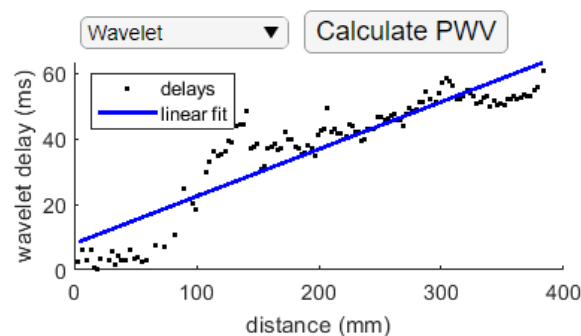

Resulting PWV can be saved by clicking Save. This records results in a matlab data struct, the PWV written to an Excel spreadsheet, and screenshot of the tool.

FUTURE to-dos:

- Read in other data formats (MATLAB .mat, Siemens 4D flow dicoms, GE PC-VIPR from UW Madison)
- Interactive angio cropping
- Refine automatic segmentation